


# ELECTRODE PLACEMENT GUIDELINES

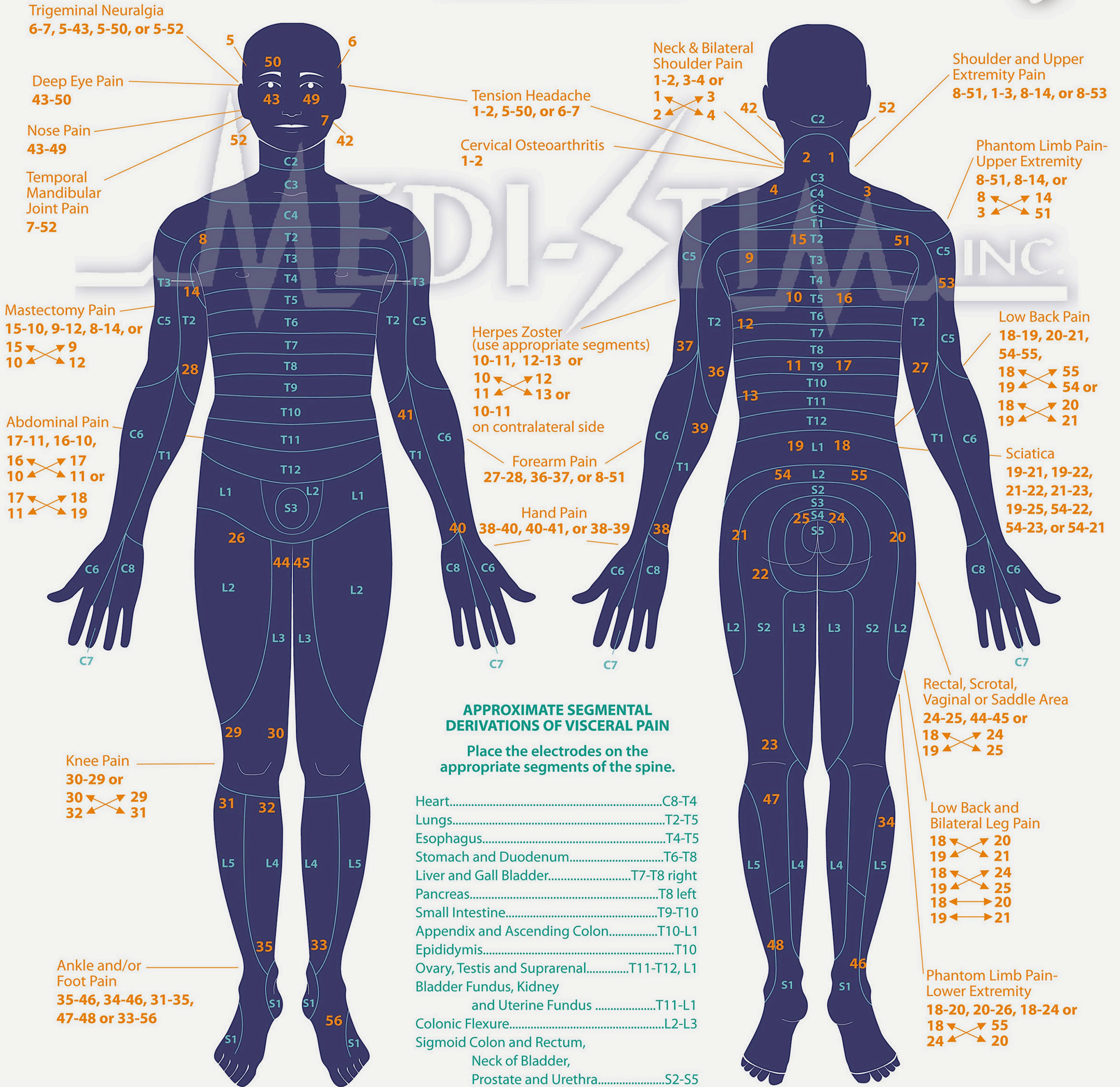

## APPROXIMATE SEGMENTAL DERIVATIONS OF VISCERAL PAIN

Place the electrodes on the appropriate segments of the spine.

Heart.....	C8-T4
Lungs.....	T2-T5
Esophagus.....	T4-T5
Stomach and Duodenum.....	T6-T8
Liver and Gall Bladder.....	T7-T8 right
Pancreas.....	T8 left
Small Intestine.....	T9-T10
Appendix and Ascending Colon.....	T10-L1
Epididymis.....	T10
Ovary, Testis and Suprarenal.....	T11-T12, L1
Bladder Fundus, Kidney and Uterine Fundus.....	T11-L1
Colonic Flexure.....	L2-L3
Sigmoid Colon and Rectum, Neck of Bladder, Prostate and Urethra.....	S2-S5